



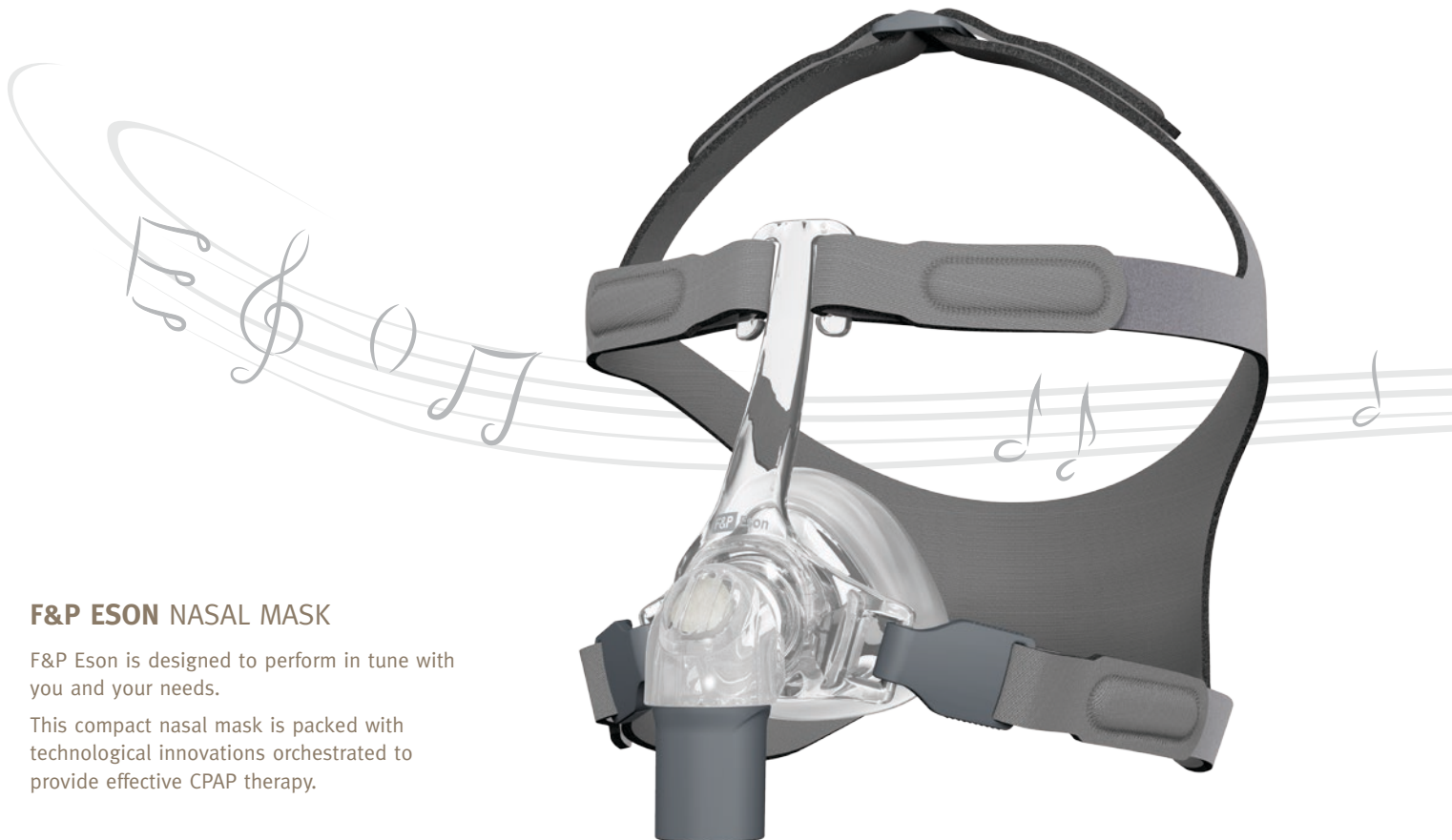

Designed to Perform – in Tune with You

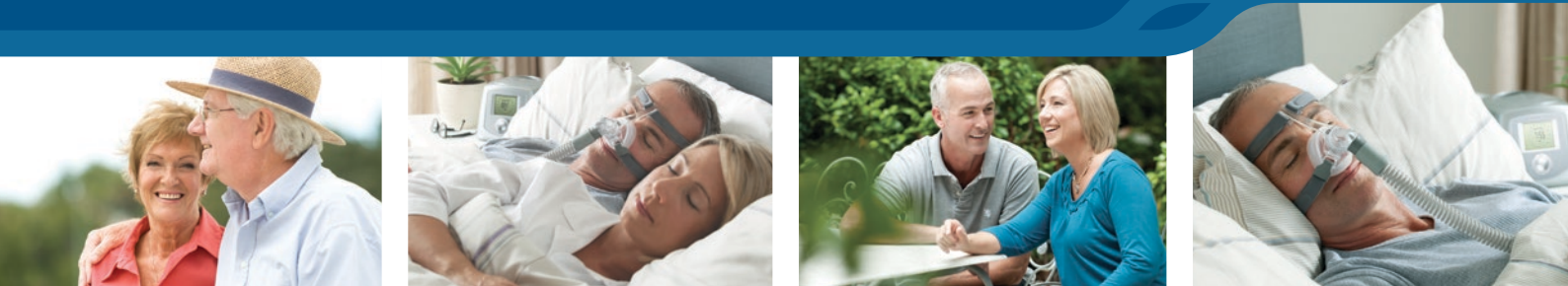


### F&P ESON NASAL MASK

F&P Eson is designed to perform in tune with you and your needs.

This compact nasal mask is packed with technological innovations orchestrated to provide effective CPAP therapy.

F&P Eson's three key components, the **RollFit™ Seal**, **ErgoFit Headgear** and **Easy Frame** work in harmony. These three key components have been designed to deliver the comfort, seal and easy use that Fisher & Paykel Healthcare masks are known for.

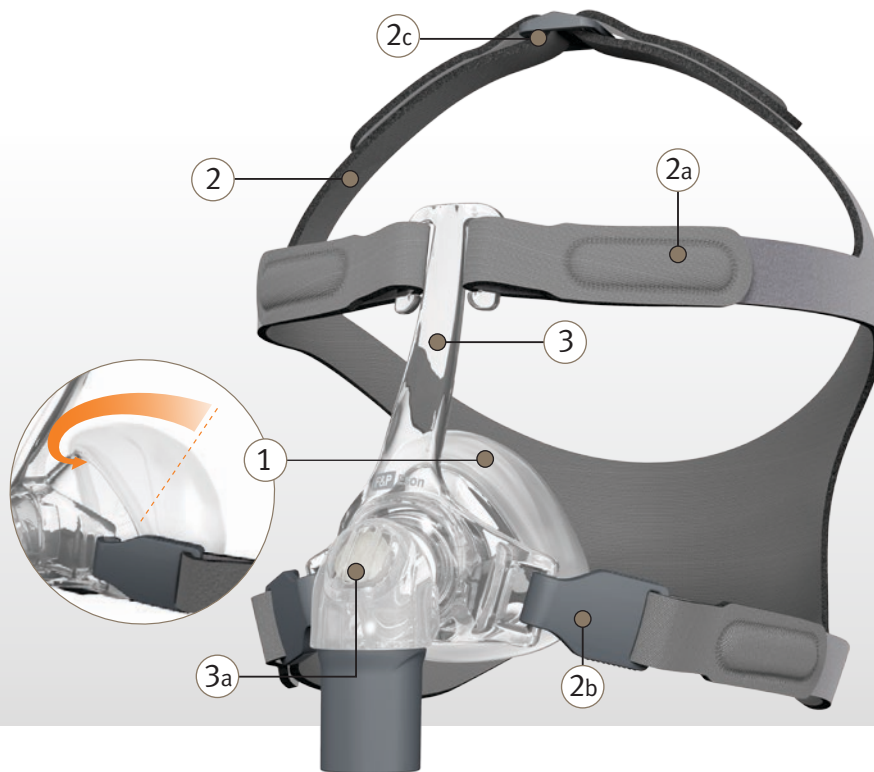

### 1. RollFit Seal

As the name suggests, the one-piece seal 'rolls' back and forth on the bridge of the nose to adjust automatically. In doing so, the **RollFit Technology** reduces pressure on the bridge of the nose.

The **RollFit Seal** incorporates **Integrated Stability Panels** to enhance sealing performance and stability.

### 2. ErgoFit Headgear

This breathable, ergonomically designed headgear has been designed to self-locate high on the rear of the head allowing for head movement (sideways, up and down).

The **ErgoFit Headgear** also incorporates embedded **VELCRO®\*** **Tabs** (2a) with **Easy-Clip hooks** (2b) for fast fitting and removal.

The adjustable **Crown Strap** (2c) enhances mask stability.

\* VELCRO® is a registered trademark of Velcro Industries B.V.

### 3. Easy Frame

This low-profile frame is stable, durable and small, and ensures a clear line of sight.

The one frame fits all three seal sizes and has an **Easy-Clip** frame attachment to aid assembly after cleaning.

The **Easy Frame** comes with the **Q Cover** (3a) for quiet use<sup>1</sup> and minimal partner draft.

The Ball and Socket Elbow rotates for freedom of movement while reducing drag from the CPAP tube.

Mask Reorder Codes		Spare Parts Reorder Codes			
400449	F&P Eson Mask Small	400HC567	F&P Eson Headgear - Small	400HC228	F&P Eson Diffuser (10pk) and Cover
400450	F&P Eson Mask Medium	400HC568	F&P Eson Headgear - Medium/Large	400HC570	F&P Eson Swivel
400451	F&P Eson Mask Large	400HC569	F&P Eson Headgear - Clips & Buckle	400HC571	F&P Eson Mask No H/G - Small
		400HC126	F&P Eson Seal - Small	400HC572	F&P Eson Mask No H/G - Medium
		400HC127	F&P Eson Seal - Medium	400HC573	F&P Eson Mask No H/G - Large
		400HC128	F&P Eson Seal - Large		